



ALIX PLIER

SIMPLIFIES A DIFFICULT TECHNIQUE FOR KW 0.8 MM TO 2 MM.

ALIX PLIER

REF.: 0227 01 00




- » Guaranteed for 5 years
- » 1 year for the jaws.
- » Stainless steel.
- » Copyright ©

Wire cutter: 1.2 mm \varnothing max

Length: 16.5 cm - Weight: 116 g

ALIX cuts by shearing (like a cigar cutter) and not by crushing. As a result, the cutting surface on the end of the pin is clean and burr-free, unlike the surface with standard cutting instruments.

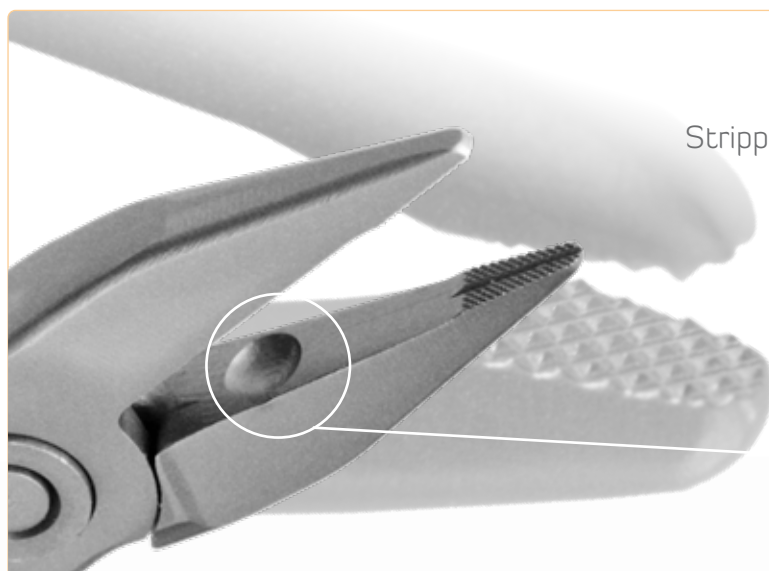

Stripped lower jaw with "V" shape groove.

To crimp KW protection tips.

DISTRIBUÉ PAR

Aked

44 C, rue de Bray
F 35510 Cesson-Sévigné

Tél. : +33 (0)2 23 35 53 43

Fax : +33 (0)9 72 26 96 49

Mail : contact@aked.bzh

AKED.bzh